


Eye Care

Model:Y401

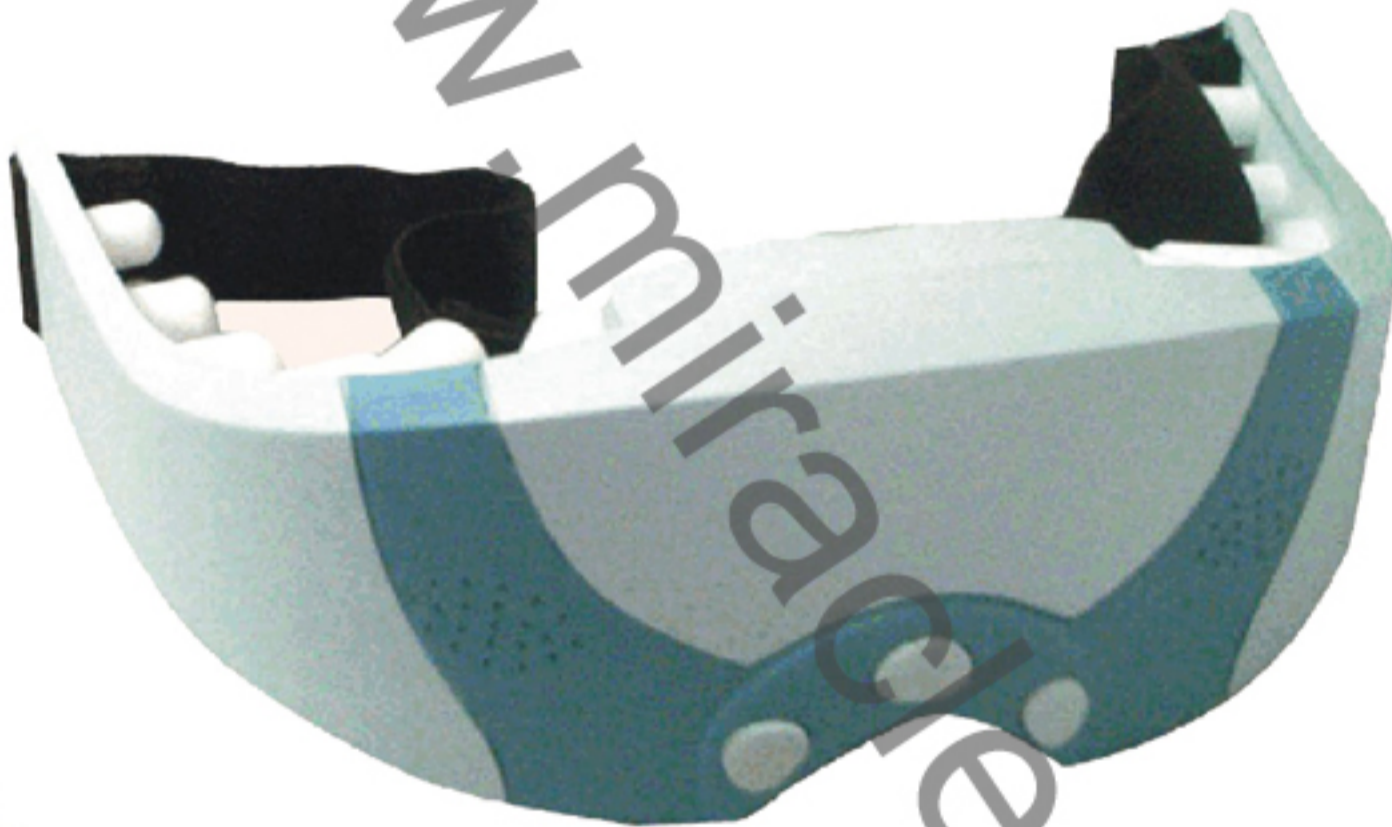

**Use it every day,
bring you bright vision**

Thanks for purchasing Eye Care Massager,
please read the instruction manual carefully before using

Using Time and Function

With the development of the borderline science subject of modern biological magnetism, plenty of research and study prove that magnetic field can have great effect on human body's tissue, organic, nervous system, biological enzyme and biological magnetic fields. Magnetic field, in the front of physical energy, when applied on special acupuncture points, can activate the function of the cell, enlarge the blood capillary and raise the level of oxygen supply, improve the nutrition state of tissue of cells, and balance the self-disciplined nerves. Therefore, it results in many kinds of medical and healthcare effects.



Acupuncture point object treatment

Zanzhu	headache, eye gall, blurred vision, walleye
Yuyao	
Sizhukong	hemicrania, remove obstruction in channels to clear eyesight
Temple	headache, eye disease
Tongjiao	refractive error, optic atrophy
Qinghou	dispel wind and remove heat, moisten and clear eyes
Chengqi	acute and chronic keratitis, lachrymation induced by wind, retinitis
Qingming	dispel wind and remove heat, moisten and clear eyes
Sibai	facial paralysis, trifacial neuralgia
Bitong	refractive error, rhinitis and the snuffle

People suitable to use the eye-massager

- * University and middle school students who suffer from false near-sightedness, eyes fatigue or eye uneasiness caused by overuse of brain or eyes through constant wearing of glasses.
- * People who work with eyes used most of the time, including computer operators, accountants, designers, secretaries, teachers, drivers driving for a long time and practice equipment users.
- * People who suffer from presbyopia and who shed tears while facing wind.
- * People who have wrinkles around the eye corners, have black eyes or stay up all night constantly.
- * People who suffer from eye-ache, headache, ache of facial nerves, neurasthenia and insomnia.

Theoretical Basis

Prevent Myopia , Beautify Eyes

Use it three times every day ,3-5 minutes every time .It enhances blood and micro-circulation around the eyes ,increase tissue oxygen uptake , thus invigorating and revitalizing the optic nerve cells ,activating key body enzymes to restoring ciliary capability ,giving you clear and sharp vision

Prevent and Cure Computer Eyes

Syndrome ,restrain neurasthenic and insomnia , remove eye fatigue

Biological magnetism science Research and Clinical studies have proven that magnetic field can have great effect on human body's tissue ,organs , nervous system ,biological enzyme and biological magnetic field ,in the form of physical energy while applied on specialised acupuncture points , can activate the tissue of cells ,enlarge the blood capillaries and raise the level of oxygen supply ,as well as improve the nutrition state of the cells , balance the self-disciplined nerves .It results in many kinds of medical and health-care effects .Applying the above theory and combining the curing theory of channel acupuncture in traditional Chinese medicine this product uses both the magnetic and medical acupuncture to activate the important acupuncture points around the eyes and harmonise the blood Qi .Thus to improve the adjusting function of the eyes muscles ,eye nerves and reach the purpose of curing near-sightedness and wrinkles and many kinds of eye-concerning diseases .;



It is proved by clinical experiment that

Eye Care Massager treatment principle illustration

When meet problem

After using :

1. Unfasten the bond .
2. Wipe the Nd.Magnetic and soft massage feeler ,keep it clean.
3. Keep it away from high temperature ,moist place.
4. Keep it away from children's touch.
5. Please take out the batteries if long time no use.

When meet problem.

Problem	Possible reason	Solution
Power off	If the batteries are used up	use new battery
	If the batteries are put in right direction	put battery correctly
Function code can't change	If the batteries are used up	use new battery
	If operate the function button exactly	re-operate it
Weak vibration	If the batteries are used up	use new battery
Can't set time	If operate the timer button exactly	re-operate it
Ache in any point	If use too long time	chose relevant use time
	If bond the elastic band too tight	adjust the elastic band
	If Nd.Magnet fall down	insert the Nd.magnet again
Soft massage feeler change color	If keep it clean and maintain it	Wipe it with wet cloth after using
	If any cosmetic left on it	Wipe it with wet cloth after using


Instructions for use

Instructions for use

- *Please take off the contact lenses before using .
- *Please keep the soft massage feeler clean (if it's dirty ,please wipe it with wet cloth).



Connect the adaptor to the socket correctly

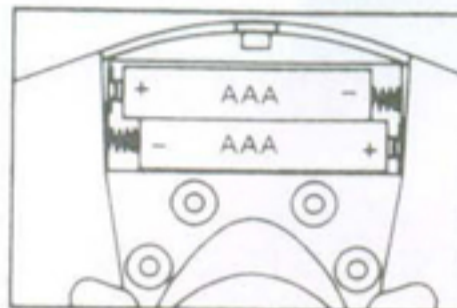

put the batteries inside the case with right direction, close the battery cover.

Instructions for use

1. Connect the adaptor or put the batteries inside ,press “on” switch, the product will work ordinally and repeatedly from Function 1to Function 9.Strong massager-Moderate Massager-Soft massager-Knead-Soft knead-Moderate knead-Strong beat-Soft beat-Moderate beat.
- 2 If press “on/off” switch again ,it will stop working.
- 3.Press “function code ” each time ,it will work according to your each choice From 1-9,they are Strong massager-Moderate massager-Soft massager-Knead-Soft knead-Moderate knead-Strong beat-Soft beat-Moderate beat.
- 4.Press “timer” ,3minutes or 5 minutes as option to set your using time.

About Nd.Magnet



If the Nd.magnet fall down from soft massage feeler, please insert it again (as attached illustration) Please switch off the product before inserting the Nd. Magnet.

Caution


Please stop using this product if you :

- *eye be injured or inflammation of eyes.
- *Head be mounded.
- *Burn around eyes or after eyes operation.

Caution:

1. Please follow the doctor's advice .
2. People who is with high level short-sight must use it under doctor's advice .
3. Children under 8 age must use it under adult's supervise.
4. Don't walk out of the door when use this product .
5. Don't use it when driving.
6. Please stop using this product after eye operation or if eye be injured.
7. Please take off the contact lenses before using .
8. Please don't use it exceeding 20minutes one time.
9. Please don't dismantle the rare-earth magnetic.
10. Please don't pluck the soft massage feeler.
11. Don't fall the product to the ground ,or immerge it in water .
12. Don't use it in bathroom or any place with high temperature.
13. It's operated by 2 "AAA" batteries or compatible adaptor .
14. Don't put the batteries with wrong direction.
15. Don't rebuild ,disconnect ,repair this product by yourself.

Part Name

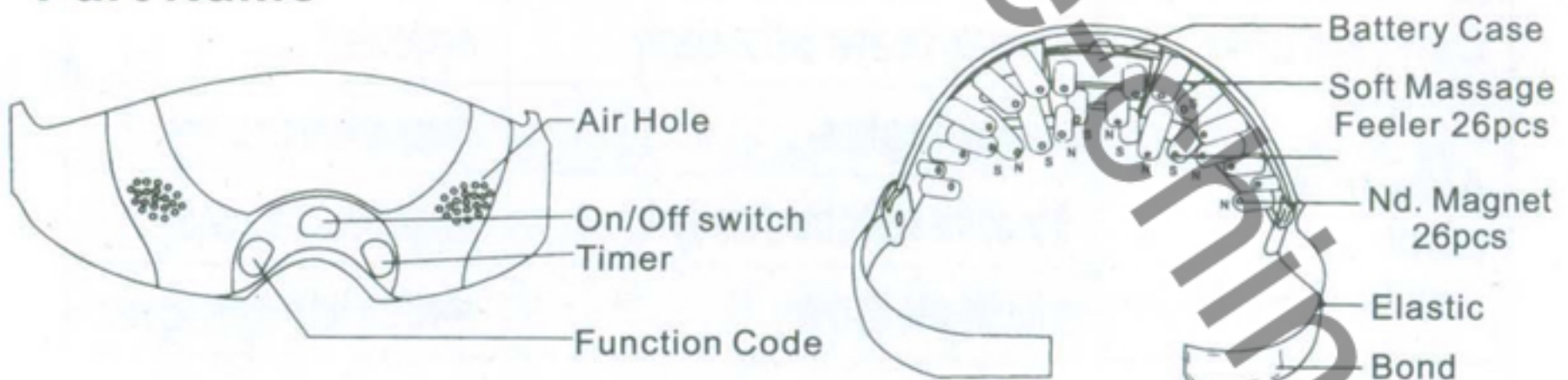

MIRACLE CO.,LTD

Address: Room B601 Hengsheng Building B , No. 205 Huanshi Zhong Road Guangzhou
China /POST NO.:510010

TEL:86-20-86673107 FAX:86-20-86672750

Email: admin@miracle-china.com

Skype: miracle-china

MSN:china.miracle@hotmail.com QQ:286372662

<http://www.miracle-china.cn>

<http://www.miracle-china.cn>